

# MediPanel

## Bed Head Panel-WP220



The unique design of modular extruded aluminium section provides a wide range of configurations to suit all needs.

A fully modular system that allows easy upgrading by adding modules to the existing panel. This minimises down time during upgrading as there is no major work interruption involved.

The design incorporates the rail system for easy support of equipment and accessories.

The wall panel system can be supplied with an optional floor rail, of fully compatible rail design.

The modular front plates, top cover and pipeline cover can be easily removed by service technicians to gain access for maintenance.

Johnson Medical's unique bed lamp has a dual movement casing which allows for it to be flipped down or flipped up. It can be used as an examination light or a night lamp.

Separate electrical and gas compartments comply with international safety standards. System is available in choice of colors to coordinate with room color.

#### Optional Accessories:

Nurse Call / TV / Telephone Modules - Additional Electric Module - Gas Module - Basket (accessory / catheter) - Cabinets - ELCB / MCB + Module - Railings - Smart Clamp™ - Room / Examination / Reading Light - Suction Cup and Accessories Tray



Bed Head Panel-WP220



Rail - Support System


Light : Flip up

Light : Flip down

*Safety through innovation™*

# JOHNSON MEDICAL

# Technical Data

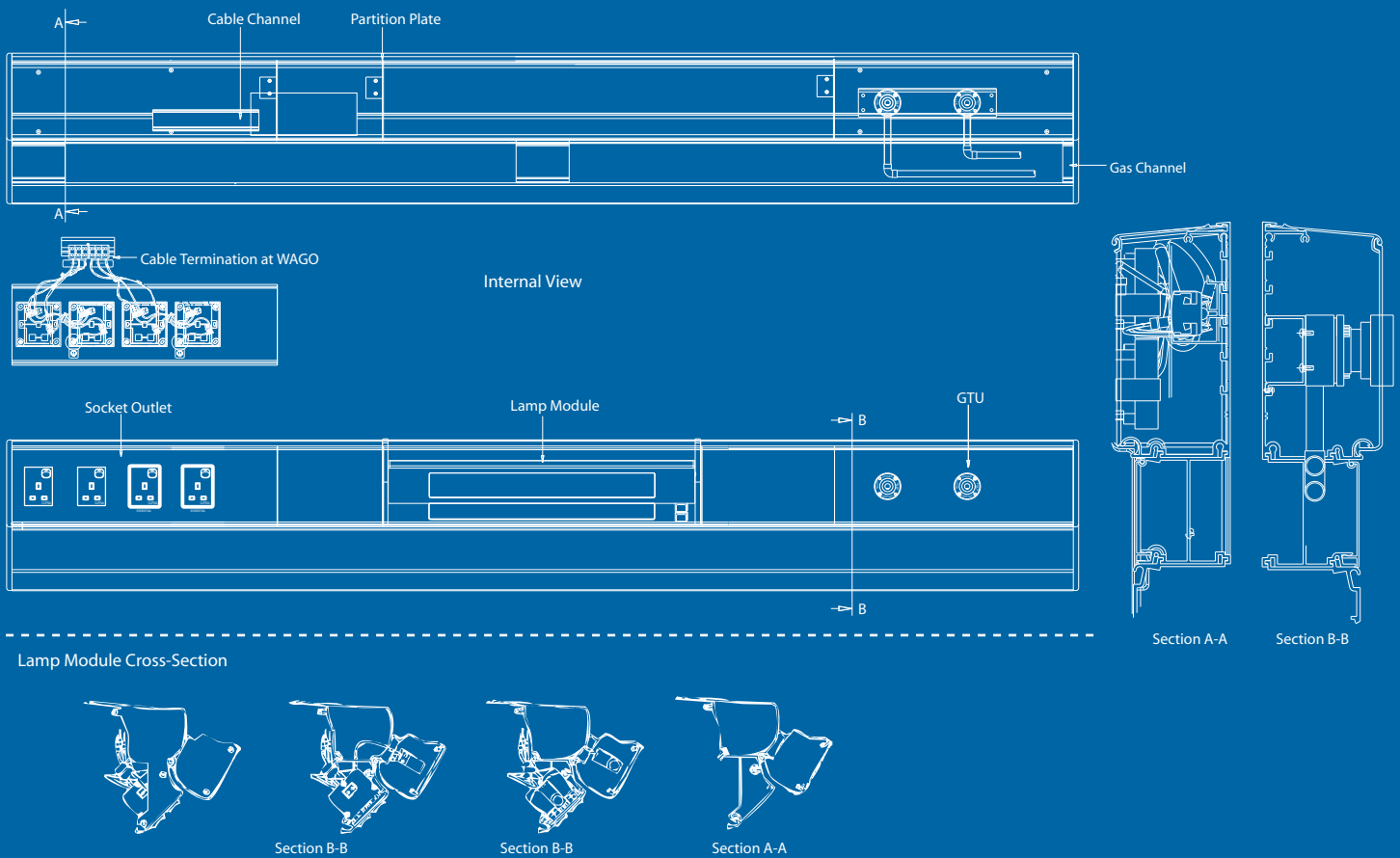

Disclaimer  
 1. These are typical configurations. They may vary from project to project based on customer requirements.  
 2. There may be improvements to the design without prior notice.

## BED HEAD PANEL - WP 220

- Major Components : Base panel and different modules made of aluminum material. Comes complete with an integrated rail system for support attachment of accessories.
- Material : Modular Aluminum sections, epoxy powder coated with a choice of colours
- Dimensions : 2000mm(l) x 285mm(w) Standard (option for customized length)
- Electrical Outlet : 120V/240V from ceiling to designated compartments. Connections are on Din-rail mounted WAGO terminal block system.
- Gas Terminal Units : Gas distributed via copper pipes mounted on the lower portion of the panel. NIST connectors used for connection
- Compliance : ISO 9001:2003, ISO 13485:2003, HTM 2022
- Product Certification : CE 120

Safely through innovation™

**JOHNSON**  
MEDICAL



HQ: Johnson Medical International Sdn Bhd

No. 12, Jalan Petaling Utama 8, 46000 Petaling Jaya, Selangor, Malaysia. Telephone: +6 03 - 7781 3081 Fax: +6 03 - 7781 9081

Email: [info@johnsonmedical.com](mailto:info@johnsonmedical.com) Website: [www.johnsonmedical.com](http://www.johnsonmedical.com)

Johnson Medical Development Pte Ltd

190, Middle Road, #12-07, Fortune Centre, 188979 Singapore. Telephone: +65 - 6338 1183 Fax: +65 - 6337 2785

Johnson Medical Representative Office

Room E, 19F, Hua Du Mansion, 838 Zhang Yang Road, Shanghai P.R. China 200122. Telephone: +86-21 5820 9723 Fax: +86-21 5081 1147